



FreeStyle  
Libre 2

## LOGIN GUIDE

This login guide provides the specific login information for your device. Please refer to the manual for all information on using myDiabby Healthcare\*.


## Step 1: Putting the FreeStyle Libre 2 data on the LibreView software

- If the monitor is scanned with the FreeStyle meter:
  - > Download the data from the meter to the LibreView software ([www.libreview.com](http://www.libreview.com)) using a computer and the yellow USB cable supplied with the meter.
- If the monitor is scanned with a smartphone and the LibreLink application, go directly to Step 2

## Step 2: Sharing LibreView data on myDiabby


From the web browser of a computer  
OR a smartphone:

- 1 Go to **www.libreview.com** and **log in** to the patient account with the **same login details** as for the LibreLink mobile app, or the LibreView login details directly.
- 2 On the LibreView website, click on:
  - The **"Curve"** symbol at the top left
  - **"Download blood sugar data"**

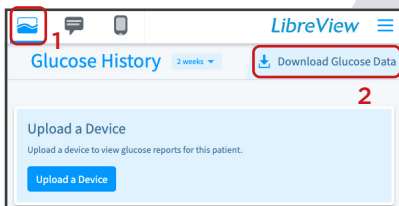


> A .csv file containing the FreeStyle Libre 2 results is saved on the phone or computer.


- 3 Open the **myDiabby** application (on smartphone OR computer) and click on **"Add data"**.

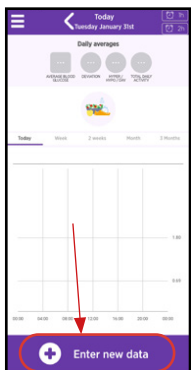


- 4 Scroll to the bottom of the screen and click on **"Import FreeStyle Libre .csv file"**.

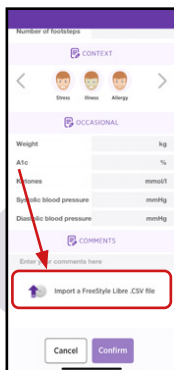

- 5 Choose the .csv file that has been saved on the phone or computer from the LibreView website

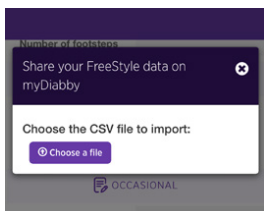

- 6 The results of the FreeStyle Libre 2 are now shared on the myDiabby app

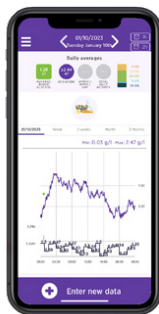

## PLEASE NOTE

All data entered on the FreeStyle Libre 2 meter and LibreLink application is shared on myDiabby when shared: Continuous glucose data, monitor scans, capillary blood glucose, insulin doses, carbohydrates, physical activity.

If you have any questions about the myDiabby application:  
[support@mydiabby.com](mailto:support@mydiabby.com) or by phone

France : 01 76 40 01 78

Belgique : 02 320 11 96



MDHC SAS  
66 avenue des Champs Elysées  
75008 Paris  
FRANCE

\*The manual is available for download and consultation on the myDiabby Healthcare platform.

