


3 Open the **myDiabby** application (on smartphone OR computer) and click on «**Enter new data**».



4 Scroll to the bottom of the screen and click on «Import a .csv file» then «Import a CareLink .csv file».

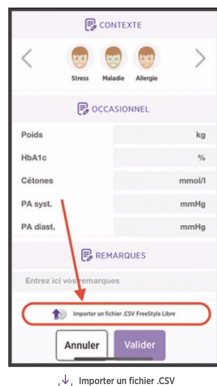

5 Select the .csv file that has been saved on your phone or computer from the CareLink site

6 **Results from the MiniMed 740G or 780G pump and its sensor are now shared on the myDiabby application.**



NOTES

All CareLink data is shared on myDiabby Healthcare when it is shared: continuous glucose data, capillary blood glucose levels, Basal, Bolus, carbohydrates, pump information, etc.

If you have any questions about the myDiabby application:
support@mydiabby.com or by phone
 France : 01 76 40 01 78 | Belgique : 02 320 11 96

*The manual is available for download and consultation on the myDiabby Healthcare platform.



Document applicable from version 2.24 of myDiabby Healthcare

myDiABBY HEALTHCARE

CONNECTION GUIDE

This connection guide provides the specific login information for your device. Please refer to the manual for all information on using myDiabby Healthcare*.




Medtronic MiniMed 780G + 740G

Phase 1 : Put pump data on the CareLink platform

If you use your pump with the **MiniMed Mobile** application, your data is automatically saved on CareLink.



If you are not using the application, download the pump data from CareLink using the reader associated with the pump. More information at: carelink.minimed.eu/app/install-uploader

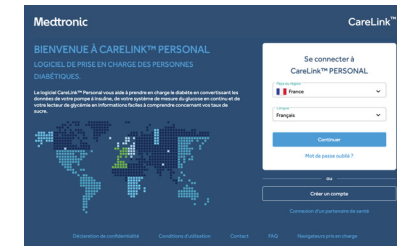


Phase 2 : Sharing CareLink data on myDiabby Healthcare

GCMCLV03D020240705EN

Method 1: Connecting your Carelink account to myDiabby Healthcare account

1 Open the **myDiabby** application (on smartphone OR computer) and click on the «**My Connected Devices**» tab.

2 Click on the «**Connect a health app**» category

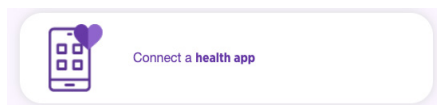

3 From the list of compatible applications, select «**CareLink**».



4 Indicate your **CareLink identifier** in the dedicated pop-up window. This is the same login you use to connect to **CareLink** or the **MiniMed Mobile** application.

Share your CareLink data on myDiabby Healthcare

Enter your CareLink username:

This is the identifier you use on CareLink or the MiniMed Mobile application

Confirm

5 Go to your **CareLink account** to accept the sharing by clicking on the link in the pop-up or via the automatic e-mail you have received.

Clinic Sharing

New requests:

myDiabby Healthcare Pending request

France

75008, FR


Patient Name: test en

Approve
Reject

Your medical provider has requested to link to your CareLink™ Personal account. When the linking process is complete, your medical provider will be able to see data from your linked medical devices that have been uploaded into your CareLink™ Personal account outside of your face to face consultations.

6 Once you've accepted the sharing, **go back to your myDiabby Healthcare account and the data will now be shared automatically!**

It may take a few minutes for the data to be recovered when you first connect.

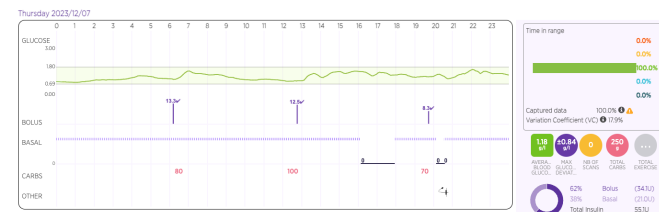

NOTES

All CareLink data is shared on myDiabby Healthcare when it is shared: continuous glucose data, capillary blood glucose levels, Basal, Bolus, carbohydrates, pump information, etc.

Method 2 : Recover data manually from CareLink

From the web browser on a computer OR a smartphone:

1 Go to carelink.minimed.com and **connect** to the patient account using the same identifiers as for the **MiniMed Mobile** application, or the **CareLink identifiers** directly.

2 On the CareLink website, click on :
- The «**Reports**» tab at the top left (curve icon)

- Then, at the bottom of the page, «**Export data (CSV)**»

> A .csv file containing the pump results is saved on the phone or computer.

